

SEPTOJECT®

SEPTOJECT NEEDLES
SEPTOJECT XL
SEPTOJECT EVOLUTION
SYRINGES



If you think a needle is just a needle, think again.

**YOUR PATIENTS FEAR NEEDLES.
SHOULDN'T YOU MINIMIZE THEIR
DISCOMFORT BY CHOOSING SEPTODONT?**



High grade surgical stainless steel tubing



Siliconized cannulas to ensure smooth insertion



Bevel indicators to help you orient the needle for different injections



Pre-threaded plastic hubs for easier threading and no metal to metal contact.

No one offers you more choices.

All our needles offer:

- Pre-threaded plastic hubs for easier threading and no metal to metal contact
- Siliconized cannulas to ensure smooth insertion
- Bevel indicators to help you orient the needle for different injections
- High grade surgical stainless steel tubing

We are the only manufacturer offering you three options:


- Septoject Needles – sharpness you can see
- Septoject XL – less pressure for the same flow
- Septoject Evolution – greater control, less discomfort

SEPTOJECT NEEDLES

TRIPLE BEVEL ►

PLASTIC HUB – EASIER THREADING

Our **Septoject Standard needles** have a triple bevel design. Many other needles on the market offer only a double bevel design. The added bevel provides for less tissue displacement and sharpness you can see.

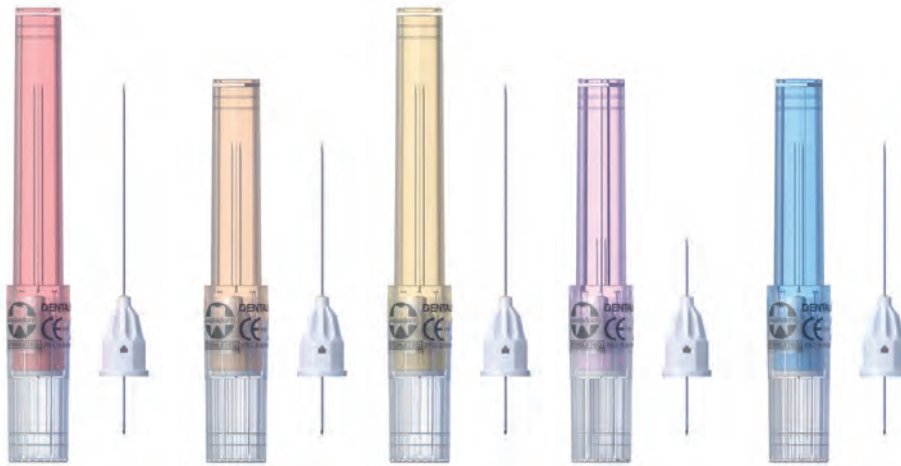

Septoject Standard

Item No.			
01N1252	25G LONG	35MM	●
01N1271	27G SHORT	25MM	●
01N1272	27G LONG	35MM	●
01N1300	30G EXTRA SHORT	10MM	●
01N1301	30G SHORT	25MM	●

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box.
10 Boxes per Case.

SEPTODONT ASPIRATING SYRINGES

Septodont dental aspirating syringes are lightweight, durable and contain no removable parts. Best of all, they are available in three distinct sizes to fit your specific injection style. They offer lightweight aluminum construction of thumb rings, wings, and a surgical grade stainless-steel body.

Septodont Standard (gold) syringe is equipped with a larger thumb ring and offers a standard length plunger for improved comfort.

Septodont Fusion (titanium) has a medium thumb ring and a standard length plunger.

Septodont Petite (blue) syringe has a compact design with a finger grip and a standard length plunger for better control and accuracy.



SEPTOJECT XL NEEDLES

LARGER LUMEN ▶

LARGER BORE – INCREASED ASPIRATION

Our **Septoject XL needles** offer a large bore design – 43% wider than a standard needle. It requires less pressure for injection.

Choose block or infiltration and simplify your selection.



Septoject XL

Item No.			
01N1500	SEPTOJECT XL - INFILTRATION	25MM	●
01N1550	SEPTOJECT XL - BLOCK	35MM	●

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box.
10 Boxes per Case.

SEPTODONT ECO SYRINGES

Our Eco syringes share the same basic design of our Septodont syringes but replace the aluminum thumb ring and wings with stainless steel making them more cost-effective.

is designed for small hands.
n provides a smaller thumb ring,
horter plunger length. This provides
added comfort when aspirating.



01-N2100

01-N2120

01-N2210

01-N2200


SEPTOJECT EVOLUTION NEEDLES

SCALPEL BEVEL ▶

GREATER CONTROL – REDUCED DISCOMFORT


Our **Septoject Evolution needle** has an innovative patented scalpel-designed bevel. It is engineered for a smoother penetration with less tissue displacement. It requires less force, provides reduced deflection and brings you better control and accuracy. Not intended for block injections.

REDUCED DEFLECTION BETTER ACCURACY

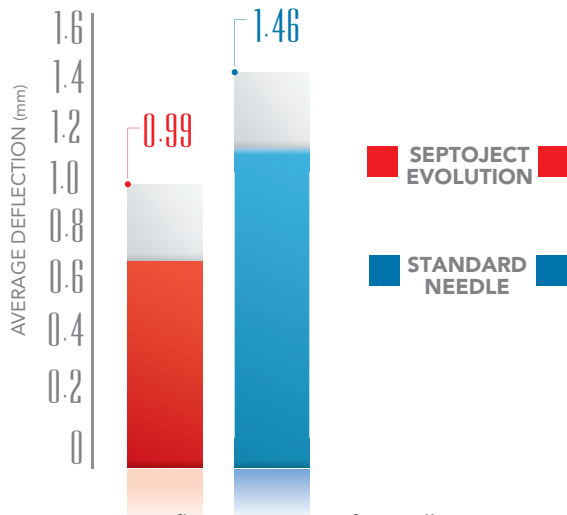

Comparitive Deflection (mm) of Needle After Penetration Through a 10mm Block of Silicone²

EASIER INSERTION - 29% LESS FORCED REQUIRED

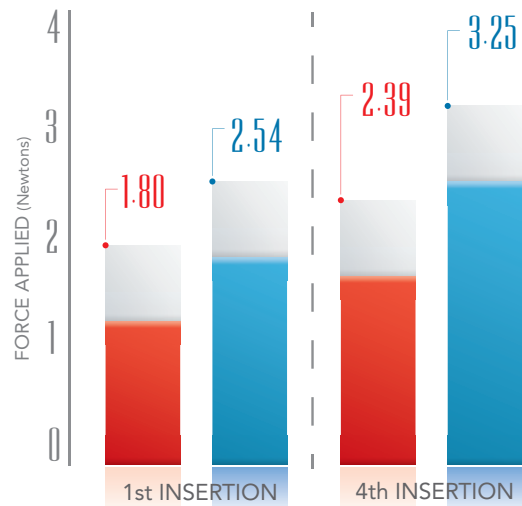

Force required (Newtons) for a 10mm Penetration of a Silicone Block¹



INFILTRATIONS

- Infiltrations (periapical)
- PSA - Posterior Superior Alveolar nerve block
- MSA - Middle Superior Alveolar nerve block
- Palatal Anesthesia



INTRALIGAMENTARY INJECTIONS

- Intraligamentary Injections
- Intraseptal Injections

Septoject Evolution

Item No.	Gauge	Length	Color
01N1600	27G SHORT	25MM	Orange
01N1630	30G SHORT	25MM	Blue

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box. 10 Boxes per Case.

Septoject Evolution

Item No.	Gauge	Length	Color
01N1650	30G EXTRA SHORT	9MM	Purple

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box. 10 Boxes per Case.



205 GRANITE RUN DRIVE
SUITE 150
LANCASTER, PA 17601

800.872.8305

www.septodontusa.com

SEPTOJECT NEEDLES
SEPTOJECT XL
SEPTOJECT EVOLUTION
SYRINGES

- | ANESTHETICS
- | PAIN CONTROL
- | ENDODONTICS
- | RESTORATIVES
- | INFECTION CONTROL
- | PERIO & SURGICAL
- | PROSTHETIC