



NEEDLES, SYRINGES

SEPTOJECT® NEEDLES

SEPTOJECT® XL

SEPTOJECT® EVOLUTION

SEPTODONT SYRINGES

From the manufacturer of Septocaine®



IF YOU THINK A NEEDLE IS JUST A NEEDLE, THINK AGAIN.

**PATIENTS FEAR NEEDLES. SHOULDN'T YOU MINIMIZE
THEIR DISCOMFORT BY CHOOSING SEPTODONT?**



Siliconized cannulas
ensure smooth insertion

High grade surgical
stainless steel tubing



NO ONE OFFERS YOU MORE CHOICES.

All our needles offer:

- Pre-threaded plastic hubs for easier threading and no metal to metal contact
- Siliconized cannulas to ensure smooth insertion
- Bevel indicators to help you orient the needle for different injections
- High grade surgical stainless steel tubing

We are the only manufacturer offering
you three options:

- **Septoject Needles** - sharpness you can see
- **Septoject XL** - less pressure for the same flow
- **Septoject Evolution** - greater control, less discomfort

Bevel indicators to help you orient
the needle for different injections

Pre-threaded plastic hubs for easier
threading and no metal to metal contact

SEPTOJECT NEEDLES // PLASTIC HUB / EASIER THREADING

Our **Septoject Standard Needles** have a triple bevel design. Many other needles on the market offer only a double bevel design. The added bevel provides for less tissue displacement and sharpness you can see.



TRIPLE
BEVEL



Item No.	Description	Length	
01N1252	25g Long	35mm	●
01N1271	27g Short	25mm	●
01N1272	27g Long	35mm	●
01N1301	30g Short	25mm	●
01N1300	30g Extra Short	10mm	●

All Gauges: 100 pre-threaded plastic hub needles per box. 10 boxes per case.

SEPTOJECT XL NEEDLES // LARGER BORE / INCREASED ASPIRATION

Our **Septoject XL Needles** offer a large bore design – 43% wider than a standard needle. They require less pressure for injection.



LARGER
LUMEN



Item No.	Description	Length	
01N1550	27g Long	35mm	●
01N1500	27g Short	25mm	●

All Gauges: 100 pre-threaded plastic hub needles per box. 10 boxes per case.

SEPTOJECT EVOLUTION NEEDLES // GREATER CONTROL / REDUCED DISCOMFORT

Our **Septoject Evolution Needles** have an innovative patented scalpel-designed bevel. Best of all, they are engineered for a smoother penetration with less tissue displacement. They require less force, provide reduced deflection and bring you better control and accuracy. Not intended for block injections.



SCALPEL BEVEL



INFILTRATIONS

- Infiltrations (periapical)
- PSA - Posterior Superior Alveolar nerve block
- MSA - Middle Superior Alveolar nerve block
- Palatal Anesthesia


INTRALIGAMENTARY INJECTIONS

- Intraligamentary Injections
- Intraseptal Injections

Item No.	Description	Length
01N1600	27g Short	25mm
01N1630	30g Short	25mm
01N1650	30g Extra Short	9mm


All Gauges: 100 pre-threaded plastic hub needles per box. 10 boxes per case.

EASIER INSERTION / 29% LESS FORCE REQUIRED

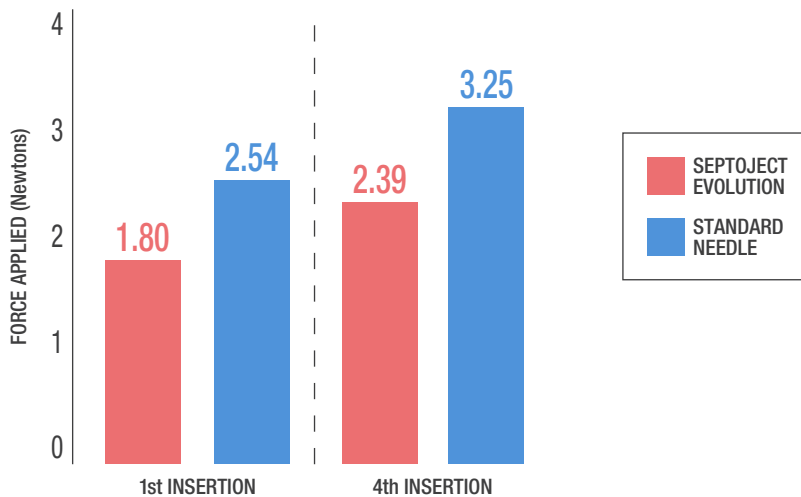

Force required (Newtons) for a 10mm penetration of a silicone block¹

REDUCED DEFLECTION / BETTER ACCURACY

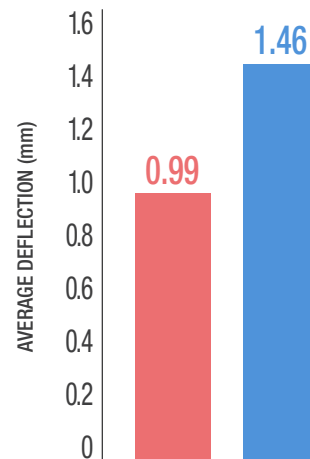

Comparative deflection (mm) of needle after penetration through a 10mm block of silicone²

1: A.Steele et al. 2: J.G. Meecan et al.

SEPTODONT SYRINGES // ASPIRATING / SELF-ASPIRATING

Septodont Dental Syringes are lightweight, durable and contain no removable parts. Best of all, they are available in three distinct sizes to fit your specific injection style. They offer lightweight aluminum construction of thumb rings, wings, and a surgical grade stainless-steel body.

Septodont Standard Syringes are equipped with a larger thumb ring and offer a standard length plunger for improved comfort. Available in **gold**.

Septodont Fusion Syringes have a medium thumb ring and a standard length plunger. Available in titanium.


Septodont Petite Syringes are designed for small hands. Its compact design provides a smaller thumb ring, finger grip and a shorter plunger length. This provides better control and added comfort when aspirating. Available in **blue** and **pink**.



Item # 01N2110
Standard Size



Item # 01N2105
Fusion Size

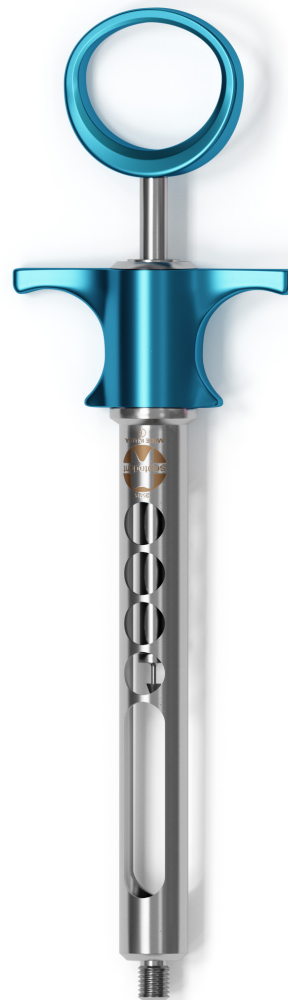


Item #01N2100
Petite Size



Item #01N2120
Petite Size

SEPTOJECT ECO SYRINGES // *COST EFFECTIVE*

Our Eco syringes share the same basic design of our Septodont syringes but replace the aluminum thumb ring and wings with stainless steel making them more cost-effective.

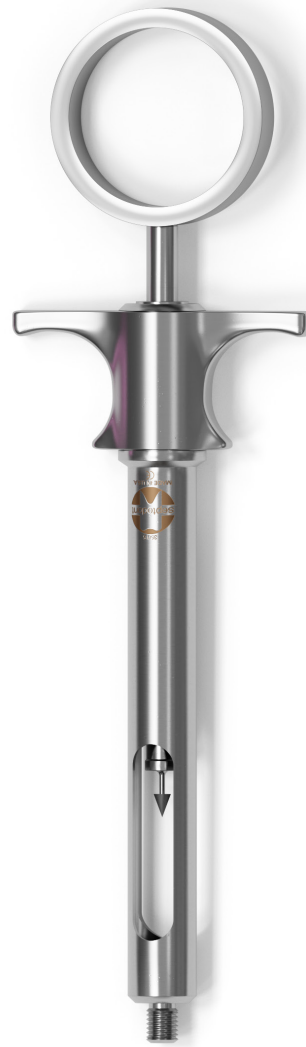

Item #01N2210
Standard Size



Item #01N2200
Standard Size

Septodont

205 Granite Run Dr., Suite 150 • Lancaster, PA 17601 • 800.872.8305

septodontusa.com

Follow us on social media channels:

R0223

